

# Schema documentation for t5\_xml.xsd

november 6, 2015

## Table of Contents

Namespace: "http://sll-mdilab.net/T5/" .....	2
Schema(s) .....	2
Main schema t5_xml.xsd .....	2
Element(s) .....	2
Element t5:PCD_01_Message .....	2
Element t5:PCD_01_Message / t5:Sending_Application .....	3
Element t5:PCD_01_Message / t5:Sending_Facility .....	4
Element t5:Patient_Result .....	5
Element t5:Patient .....	5
Element t5:Identifier .....	6
Element t5:Patient / t5:Name .....	6
Element t5:Location .....	7
Element t5:Order_Observations .....	7
Element t5:Order .....	8
Element t5:Order / t5:FillerOrderNumber .....	8
Element t5:Order / t5:UniversalServiceID .....	9
Element t5:MDS .....	10
Element t5:VMD .....	11
Element t5:CHAN .....	11
Element t5:Metric .....	12
Element t5:Metric / t5:Facet .....	13
Element t5:Observation .....	14
Element t5:Observation / t5:ObsIdentifier .....	16
Element t5:Value .....	18
Element t5:Unit .....	18
Element t5:Timestamp .....	19
Element t5:Method .....	19
Element t5:Body_Site .....	19
Element t5:Status .....	19
Element t5:Abnormal_Flag .....	20
Element t5:EquipmentIdentifier .....	20
Element t5:Note .....	20
Element t5:Facet .....	21
Element t5:Result .....	21
Complex Type(s) .....	22
Complex Type t5:PCD_01_Message .....	22
Complex Type t5:Patient_Result .....	23
Complex Type t5:Patient .....	24
Complex Type t5:Location .....	24
Complex Type t5:Order_Observations .....	25
Complex Type t5:Order .....	25
Complex Type t5:MDS .....	26
Complex Type t5:VMD .....	27
Complex Type t5:CHAN .....	28
Complex Type t5:Metric .....	29
Complex Type t5:Observation .....	30
Complex Type t5:Unit .....	33
Complex Type t5:Timestamp .....	34
Complex Type t5:EquipmentIdentifier .....	34
Complex Type t5:Facet .....	34
Namespace: "" .....	35
Attribute(s) .....	35
Attribute t5:PCD_01_Message / t5:Sending_Application / @idLocal .....	35
Attribute t5:PCD_01_Message / t5:Sending_Application / @idUniversal .....	35
Attribute t5:PCD_01_Message / t5:Sending_Application / @idUniversalType .....	36
Attribute t5:PCD_01_Message / t5:Sending_Facility / @idLocal .....	36
Attribute t5:PCD_01_Message / t5:Sending_Facility / @idUniversal .....	36
Attribute t5:PCD_01_Message / t5:Sending_Facility / @idUniversalType .....	36
Attribute t5:Identifier / @authority .....	36
Attribute t5:Identifier / @typeCode .....	37
Attribute t5:Location / @room .....	37
Attribute t5:Location / @bed .....	37

Attribute t5:Location / @pointOfCare .....	37
Attribute t5:Order / t5:FillerOrderNumber / @idLocal .....	37
Attribute t5:Order / t5:FillerOrderNumber / @nameSpaceID .....	38
Attribute t5:Order / t5:FillerOrderNumber / @idUniversal .....	38
Attribute t5:Order / t5:FillerOrderNumber / @idUniversalType .....	38
Attribute t5:Order / t5:UniversalServiceID / @id .....	38
Attribute t5:Order / t5:UniversalServiceID / @codingSystemName .....	38
Attribute t5:Order / @timeStamp .....	39
Attribute t5:Observation / t5:ObsIdentifier / @id .....	39
Attribute t5:Observation / t5:ObsIdentifier / @codingSystemName .....	39
Attribute t5:Observation / t5:ObsIdentifier / @isAlternate .....	39
Attribute t5:Observation / t5:ObsIdentifier / @notFound .....	39
Attribute t5:Value / @typeHL7v2 .....	40
Attribute t5:Unit / @id .....	40
Attribute t5:Unit / @codingSystemName .....	40
Attribute t5:EquipmentIdentifier / @nameSpaceID .....	40
Attribute t5:Observation / @uid .....	40
Attribute t5:Observation / @hierarchy .....	41
Attribute t5:Observation / @index .....	41
Attribute t5:Observation / @setid .....	41
Attribute t5:Observation / @violatesIHE .....	41
Attribute t5:Metric / t5:Facet / @index .....	41
Attribute t5:Metric / t5:Facet / @id .....	42
Attribute t5:Metric / @index .....	42
Attribute t5:Metric / @id .....	42
Attribute t5:CHAN / @index .....	42
Attribute t5:CHAN / @id .....	42
Attribute t5:VMD / @index .....	43
Attribute t5:VMD / @id .....	43
Attribute t5:MDS / @index .....	43
Attribute t5:MDS / @id .....	43
Attribute t5:PCD_01_Message / @timeStamp .....	43
Attribute t5:PCD_01_Message / @id .....	44
Attribute t5:Facet / @index .....	44
Attribute t5:Facet / @id .....	44
Attribute t5:Result / @unit .....	44

**Namespace: "http://sll-mdilab.net/T5/"**

## Schema(s)

**Main schema t5\_xml.xsd**

Namespace	http://sll-mdilab.net/T5/
Annotations	T5 observations
Properties	attribute form default: unqualified
	element form default: qualified

## Element(s)

**Element t5:PCD\_01\_Message**

Namespace	http://sll-mdilab.net/T5/
-----------	---------------------------



Properties	content: complex			
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>	
	<b>idLocal</b>		optional	
		MSH-3.1		
	<b>idUniversal</b>		optional	
		MSH-3.2		
	<b>idUniversalType</b>		optional	
		MSH-3.3		
Source	<pre>&lt;xs:element name="Sending_Application"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-3 MSH-3.1 is used as Application ID, which is an identifier of the System&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:complexType&gt;     &lt;xs:attribute name="idLocal" use="optional"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;MSH-3.1&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="idUniversal"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;MSH-3.2&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="idUniversalType"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;MSH-3.3&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt;</pre>			

### Element t5:PCD\_01\_Message / t5:Sending\_Facility

Namespace	http://sll-mdilab.net/T5/																															
Annotations	MSH-4 to distinguish Hospital/department																															
Diagram																																
Properties	content: complex																															
Attributes	<table><thead><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr></thead><tbody><tr><td>idLocal</td><td></td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">MSH-4.1</td></tr><tr><td>idUniversal</td><td></td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">MSH-4.2</td></tr><tr><td>idUniversalType</td><td></td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">MSH-4.3</td></tr></tbody></table>	QName	Type	Use		idLocal		optional			MSH-4.1			idUniversal		optional			MSH-4.2			idUniversalType		optional			MSH-4.3					
QName	Type	Use																														
idLocal		optional																														
	MSH-4.1																															
idUniversal		optional																														
	MSH-4.2																															
idUniversalType		optional																														
	MSH-4.3																															
Source	<pre>&lt;xs:element name="Sending_Facility"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-4&lt;/xs:documentation&gt;     &lt;xs:documentation&gt;to distinguish Hospital/department&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:complexType&gt;</pre>																															

```

<xs:attribute name="idLocal" use="optional">
  <xs:annotation>
    <xs:documentation>MSH-4.1</xs:documentation>
  </xs:annotation>
</xs:attribute>
<xs:attribute name="idUniversal">
  <xs:annotation>
    <xs:documentation>MSH-4.2</xs:documentation>
  </xs:annotation>
</xs:attribute>
<xs:attribute name="idUniversalType">
  <xs:annotation>
    <xs:documentation>MSH-4.3</xs:documentation>
  </xs:annotation>
</xs:attribute>
</xs:complexType>
</xs:element>

```

## Element t5:Patient\_Result

Namespace	http://sll-mdilab.net/T5/
Diagram	
Type	t5:Patient_Result
Properties	content: complex
Used by	Complex Type t5:PCD_01_Message
Model	t5:Patient , t5:Order_Observations+
Children	t5:Order_Observations, t5:Patient
Instance	<pre> &lt;t5:Patient_Result xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Patient&gt;{1,1}&lt;/t5:Patient&gt;   &lt;t5:Order_Observations&gt;{1,unbounded}&lt;/t5:Order_Observations&gt; &lt;/t5:Patient_Result&gt; </pre>
Source	<pre> &lt;xs:element name="Patient_Result" type="t5:Patient_Result"&gt; &lt;/xs:element&gt; </pre>

## Element t5:Patient

Namespace	http://sll-mdilab.net/T5/
Diagram	
Type	t5:Patient
Properties	content: complex
Used by	Complex Type t5:Patient_Result
Model	t5:Identifier{1,6} , t5:Name , t5:Location{0,1}

Children	t5:Identifier, t5:Location, t5:Name
Instance	<pre>&lt;t5:Patient xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Identifier authority="" typeCode=""&gt;{1,6}&lt;/t5:Identifier&gt;   &lt;t5:Name&gt;{1,1}&lt;/t5:Name&gt;   &lt;t5:Location bed="" pointOfCare="" room=""&gt;{0,1}&lt;/t5:Location&gt; &lt;/t5:Patient&gt;</pre>
Source	<pre>&lt;xs:element name="Patient" type="t5:Patient"&gt; &lt;/xs:element&gt;</pre>

## Element t5:Identifier

Namespace	http://sll-mdilab.net/T5/		
Annotations	PID-3.1		
Diagram			
Type	extension of xs:string		
Properties	content:	complex	
	mixed:	true	
Used by	Complex Type	t5:Patient	
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>
	<b>authority</b>		required
		PID-3.4	
	<b>typeCode</b>		optional
		PID-3.5	
Source	<pre>&lt;xs:element name="Identifier"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PID-3.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:complexType mixed="true"&gt;     &lt;xs:simpleContent&gt;       &lt;xs:extension base="xs:string"&gt;         &lt;xs:attribute name="authority" use="required"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;PID-3.4&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;         &lt;xs:attribute name="typeCode"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;PID-3.5&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;       &lt;/xs:extension&gt;     &lt;/xs:simpleContent&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt;</pre>		

## Element t5:Patient / t5:Name

Namespace	http://sll-mdilab.net/T5/
Annotations	PID-5
Diagram	

Source	<pre>&lt;xs:element name="Name"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PID-5&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:element&gt;</pre>
--------	---

## Element t5:Location

Namespace	http://sll-mdilab.net/T5/																															
Diagram																																
Type	t5:Location																															
Properties	content:	complex																														
Used by	Complex Type	t5:Patient																														
Attributes	<table><thead><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr></thead><tbody><tr><td><b>bed</b></td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>PV1-3.3</td><td></td><td></td></tr><tr><td><b>pointOfCare</b></td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>PV1-3.1</td><td></td><td></td></tr><tr><td><b>room</b></td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>PV1-3.2</td><td></td><td></td></tr></tbody></table>	QName	Type	Use		<b>bed</b>	xs:string	optional			PV1-3.3			<b>pointOfCare</b>	xs:string	optional			PV1-3.1			<b>room</b>	xs:string	optional			PV1-3.2					
QName	Type	Use																														
<b>bed</b>	xs:string	optional																														
	PV1-3.3																															
<b>pointOfCare</b>	xs:string	optional																														
	PV1-3.1																															
<b>room</b>	xs:string	optional																														
	PV1-3.2																															
Source	<xs:element name="Location" type="t5:Location" />																															

## Element t5:Order\_Observations

Namespace	http://sll-mdilab.net/T5/		
Diagram			
Type	t5:Order_Observations		
Properties	content:	complex	
Used by	Complex Type	t5:Patient_Result	

Model	t5:Order{0,1} , t5:MDS* , t5:Note*
Children	t5:MDS, t5:Note, t5:Order
Instance	<pre>&lt;t5:Order_Observations xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Order timeStamp=""&gt;{0,1}&lt;/t5:Order&gt;   &lt;t5:MDS id="" index=""&gt;{0,unbounded}&lt;/t5:MDS&gt;   &lt;t5:Note&gt;{0,unbounded}&lt;/t5:Note&gt; &lt;/t5:Order_Observations&gt;</pre>
Source	<pre>&lt;xs:element name="Order_Observations" type="t5:Order_Observations"&gt; &lt;/xs:element&gt;</pre>

## Element t5:Order

Namespace	http://sll-mdilab.net/T5/			
Diagram				
Type	t5:Order			
Properties	content:	complex		
Used by	Complex Type	t5:Order_Observations		
Model	t5:FillerOrderNumber , t5:UniversalServiceID			
Children	t5:FillerOrderNumber, t5:UniversalServiceID			
Instance	<pre>&lt;t5:Order timeStamp="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:FillerOrderNumber idLocal="" idUniversal="" idUniversalType="" namespaceID=""&gt;{1,1}&lt;/t5:FillerOrderNumber&gt;   &lt;t5:UniversalServiceID codingSystemName="" id=""&gt;{1,1}&lt;/t5:UniversalServiceID&gt; &lt;/t5:Order&gt;</pre>			
Attributes	QName	Type	Use	
	timeStamp	xs:dateTime	optional	
		OBR-7		
Source	<pre>&lt;xs:element name="Order" type="t5:Order"&gt; &lt;/xs:element&gt;</pre>			

## Element t5:Order / t5:FillerOrderNumber

Namespace	http://sll-mdilab.net/T5/
Annotations	OBR-3.1



Diagram																																					
Properties	<table><tr><td>content:</td><td>complex</td></tr><tr><td>maxOccurs:</td><td>1</td></tr></table>	content:	complex	maxOccurs:	1																																
content:	complex																																				
maxOccurs:	1																																				
Attributes	<table><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr><tr><td>idLocal</td><td></td><td>optional</td><td></td></tr><tr><td></td><td>OBR-3.1</td><td></td><td></td></tr><tr><td>idUniversal</td><td></td><td>optional</td><td></td></tr><tr><td></td><td>OBR-3.3</td><td></td><td></td></tr><tr><td>idUniversalType</td><td></td><td>optional</td><td></td></tr><tr><td></td><td>MSH-3.4</td><td></td><td></td></tr><tr><td>nameSpaceID</td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>OBR-3.2</td><td></td><td></td></tr></table>	QName	Type	Use		idLocal		optional			OBR-3.1			idUniversal		optional			OBR-3.3			idUniversalType		optional			MSH-3.4			nameSpaceID	xs:string	optional			OBR-3.2		
QName	Type	Use																																			
idLocal		optional																																			
	OBR-3.1																																				
idUniversal		optional																																			
	OBR-3.3																																				
idUniversalType		optional																																			
	MSH-3.4																																				
nameSpaceID	xs:string	optional																																			
	OBR-3.2																																				
Source	<pre>&lt;xs:element name="FillerOrderNumber" maxOccurs="1"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBR-3.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:complexType&gt;     &lt;xs:attribute name="idLocal" use="optional"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-3.1&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="nameSpaceID" type="xs:string"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-3.2&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="idUniversal"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-3.3&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="idUniversalType"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;MSH-3.4&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt;</pre>																																				

**Element t5:Order / t5:UniversalServiceID**

Namespace	http://sll-mdilab.net/T5/
Annotations	OBR-4.2 text

Diagram	<div><div>UniversalServiceID</div><div>OBR-4.2 text</div><div><div>@ Attributes</div><div><div>@ id</div><div>OBR-4.1 id</div></div><div><div>@ codingSystemName</div><div>Type xs:string</div><div>OBR-4.3</div></div></div></div>																							
Properties	content:	complex																						
	minOccurs:	1																						
	maxOccurs:	1																						
Attributes	<table><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr><tr><td><b>codingSystemName</b></td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>OBR-4.3</td><td></td><td></td></tr><tr><td><b>id</b></td><td></td><td>optional</td><td></td></tr><tr><td></td><td>OBR-4.1 id</td><td></td><td></td></tr></table>	QName	Type	Use		<b>codingSystemName</b>	xs:string	optional			OBR-4.3			<b>id</b>		optional			OBR-4.1 id					
QName	Type	Use																						
<b>codingSystemName</b>	xs:string	optional																						
	OBR-4.3																							
<b>id</b>		optional																						
	OBR-4.1 id																							
Source	<pre>&lt;xs:element name="UniversalServiceID" maxOccurs="1" minOccurs="1"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBR-4.2 text&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:complexType&gt;     &lt;xs:attribute name="id" use="optional"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-4.1 id&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="codingSystemName" type="xs:string" use="optional"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-4.3&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt;</pre>																							

## Element t5:MDS

Namespace	http://sll-mdilab.net/T5/		
Diagram			
Type	t5:MDS		
Properties	content:	complex	
Used by	Complex Type	t5:Order_Observations	
Model	t5:VMD{0,1} , t5:Observation*		

Children	t5:Observation, t5:VMD			
Instance	<pre>&lt;t5:MDS id="" index="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:VMD id="" index=""&gt;{0,1}&lt;/t5:VMD&gt;   &lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE=""&gt;{0,unbounded}&lt;/ t5:Observation&gt; &lt;/t5:MDS&gt;</pre>			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.0.0.0		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:element name="MDS" type="t5:MDS"&gt; &lt;/xs:element&gt;</pre>			

## Element t5:VMD

Namespace	http://sll-mdilab.net/T5/			
Diagram				
Type	t5:VMD			
Properties	content:	complex		
Used by	Complex Type	t5:MDS		
Model	t5:CHAN{0,1} , t5:Observation*			
Children	t5:CHAN, t5:Observation			
Instance	<pre>&lt;t5:VMD id="" index="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:CHAN id="" index=""&gt;{0,1}&lt;/t5:CHAN&gt;   &lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE=""&gt;{0,unbounded}&lt;/ t5:Observation&gt; &lt;/t5:VMD&gt;</pre>			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.y.0.0		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:element name="VMD" type="t5:VMD"&gt; &lt;/xs:element&gt;</pre>			

## Element t5:CHAN

Namespace	http://sll-mdilab.net/T5/
-----------	---------------------------

Diagram				
Type	t5:CHAN			
Properties	content:	complex		
Used by	Complex Type	t5:VMD		
Model	t5:Metric{0,1} , t5:Observation*			
Children	t5:Metric, t5:Observation			
Instance	<pre>&lt;t5:CHAN id="" index="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Metric id="" index=""&gt;{0,1}&lt;/t5:Metric&gt;   &lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE=""&gt;{0,unbounded}&lt;/ t5:Observation&gt; &lt;/t5:CHAN&gt;</pre>			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.y.z.0		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:element name="CHAN" type="t5:CHAN"&gt; &lt;/xs:element&gt;</pre>			

## Element t5:Metric

Namespace	http://sll-mdilab.net/T5/			
Diagram				

Type	t5:Metric			
Properties	content:	complex		
Used by	Complex Type	t5:CHAN		
Model	t5:Facet{0,1} , t5:Observation*			
Children	t5:Facet, t5:Observation			
Instance	<pre>&lt;t5:Metric id="" index="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Facet id="" index=""&gt;{0,1}&lt;/t5:Facet&gt;   &lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE=""&gt;{0,unbounded}&lt;/ t5:Observation&gt; &lt;/t5:Metric&gt;</pre>			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.y.z.w		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:element name="Metric" type="t5:Metric"&gt; &lt;/xs:element&gt;</pre>			

## Element t5:Metric / t5:Facet

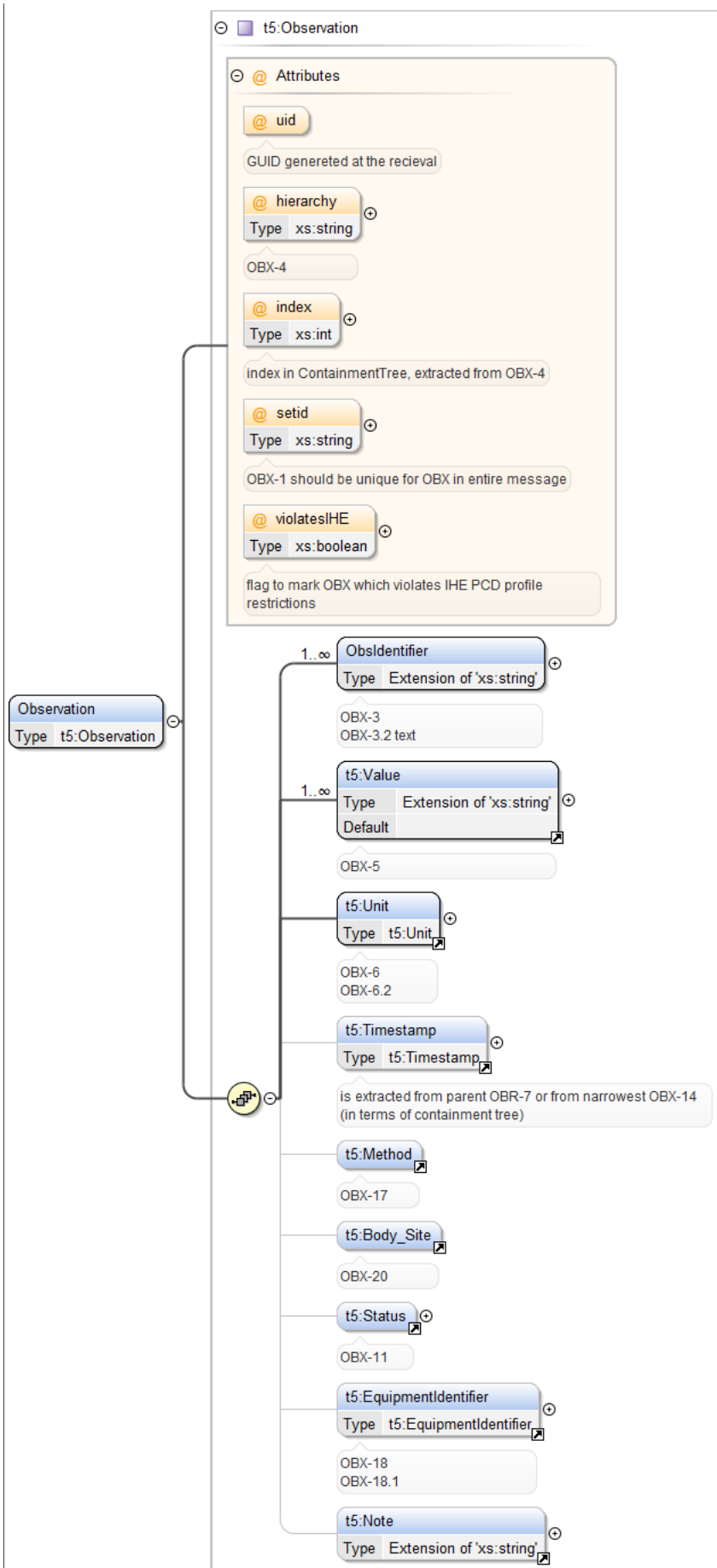
Namespace	http://sll-mdilab.net/T5/							
Diagram								
Properties	<table><tr><td>content:</td><td>complex</td></tr><tr><td>minOccurs:</td><td>0</td></tr></table>				content:	complex	minOccurs:	0
content:	complex							
minOccurs:	0							
Model	t5:Observation							
Children	t5:Observation							
Instance	<pre>&lt;t5:Facet id="" index="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE=""&gt;{1,1}&lt;/t5:Observation&gt; &lt;/t5:Facet&gt;</pre>							
Attributes	QName	Type	Use					
	id		optional					
		OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t						
	index	xs:int	optional					
		index in ContainmentTree, extracted from OBX-4						
Source	<pre>&lt;xs:element name="Facet" minOccurs="0"&gt;   &lt;xs:complexType&gt;     &lt;xs:sequence&gt;       &lt;xs:element ref="t5:Observation"/&gt;     &lt;/xs:sequence&gt;     &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:attribute&gt;     &lt;xs:attribute name="id"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t&lt;/ xs:documentation&gt;</pre>							

	<pre>&lt;/xs:annotation&gt; &lt;/xs:attribute&gt; &lt;/xs:complexType&gt; &lt;/xs:element&gt;</pre>
--	---

**Element t5:Observation**

Namespace	<a href="http://sll-mdilab.net/T5/">http://sll-mdilab.net/T5/</a>
-----------	---

## Diagram



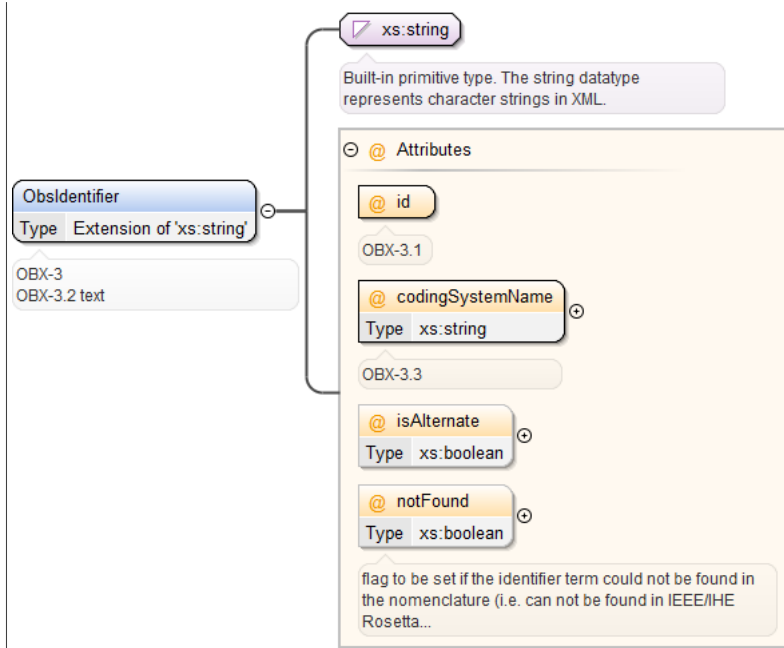
Type	t5:Observation			
Properties	content:	complex		
Used by	Element	t5:Metric/t5:Facet		
	Complex Types	t5:CHAN, t5:Facet, t5:MDS, t5:Metric, t5:VMD		
Model	t5:ObsIdentifier+ , t5:Value+ , t5:Unit , t5:Timestamp{0,1} , t5:Method{0,1} , t5:Body_Site{0,1} , t5:Status{0,1} , t5:EquipmentIdentifier{0,1} , t5:Note{0,1}			
Children	t5:Body_Site, t5:EquipmentIdentifier, t5:Method, t5:Note, t5:ObsIdentifier, t5:Status, t5:Timestamp, t5:Unit, t5:Value			
Instance	<pre>&lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:ObsIdentifier codingSystemName="" id="" isAlternate="" notFound=""&gt;{1,unbounded}&lt;/t5:ObsIdentifier&gt;   &lt;t5:Value typeHL7v2=""&gt;{1,unbounded}&lt;/t5:Value&gt;   &lt;t5:Unit codingSystemName="" id=""&gt;{1,1}&lt;/t5:Unit&gt;   &lt;t5:Timestamp&gt;{0,1}&lt;/t5:Timestamp&gt;   &lt;t5:Method&gt;{0,1}&lt;/t5:Method&gt;   &lt;t5:Body_Site&gt;{0,1}&lt;/t5:Body_Site&gt;   &lt;t5:Status&gt;{0,1}&lt;/t5:Status&gt;   &lt;t5:EquipmentIdentifier namespaceID=""&gt;{0,1}&lt;/t5:EquipmentIdentifier&gt;   &lt;t5:Note&gt;{0,1}&lt;/t5:Note&gt; &lt;/t5:Observation&gt;</pre>			
Attributes	QName	Type	Use	
	hierarchy	xs:string	optional	
		OBX-4		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
	setid	xs:string	optional	
		OBX-1 should be unique for OBX in entire message		
	uid		optional	
		GUID generated at the recieval		
	violatesIHE	xs:boolean	optional	
		flag to mark OBX which violates IHE PCD profile restrictions		
Source	<pre>&lt;xs:element name="Observation" type="t5:Observation"&gt; &lt;/xs:element&gt;</pre>			

## Element t5:Observation / t5:ObsIdentifier

Namespace	http://sll-mdilab.net/T5/
Annotations	OBX-3
	OBX-3.2 text



## Diagram



Type extension of xs:string

Properties content: complex  
maxOccurs: unbounded

Attributes	QName	Type	Use	
	<b>codingSystemName</b>	xs:string	required	
		OBX-3.3		
	<b>id</b>		required	
		OBX-3.1		
	<b>isAlternate</b>	xs:boolean	optional	
	<b>notFound</b>	xs:boolean	optional	
		flag to be set if the identifier term could not be found in the nomenclature (i.e. can not be found in IEEE/IHE Rosetta table)		

Source

```
<xs:element name="ObsIdentifier" maxOccurs="unbounded">
  <xs:annotation>
    <xs:documentation>OBX-3</xs:documentation>
    <xs:documentation>OBX-3.2 text</xs:documentation>
  </xs:annotation>
  <xs:complexType>
    <xs:simpleContent>
      <xs:extension base="xs:string">
        <xs:attribute name="id" use="required">
          <xs:annotation>
            <xs:documentation>OBX-3.1</xs:documentation>
          </xs:annotation>
        </xs:attribute>
        <xs:attribute name="codingSystemName" type="xs:string" use="required">
          <xs:annotation>
            <xs:documentation>OBX-3.3</xs:documentation>
          </xs:annotation>
        </xs:attribute>
        <xs:attribute name="isAlternate" type="xs:boolean"/>
        <xs:attribute name="notFound" type="xs:boolean">
          <xs:annotation>
            <xs:documentation>flag to be set if the identifier term could not be found in the
nomenclature (i.e. can not be found in IEEE/IHE Rosetta table)</xs:documentation>
          </xs:annotation>
        </xs:attribute>
      </xs:extension>
    </xs:simpleContent>
  </xs:complexType>
</xs:element>
```

## Element t5:Value

Namespace	http://sll-mdilab.net/T5/		
Diagram			
Type	extension of xs:string		
Properties	content:	complex	
	default:		
	mixed:	true	
Used by	Complex Type	t5:Observation	
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>
	typeHL7v2	xs:string	optional
Source	<pre> &lt;xs:element name="Value" default=""&gt;   &lt;xs:complexType mixed="true"&gt;     &lt;xs:simpleContent&gt;       &lt;xs:extension base="xs:string"&gt;         &lt;xs:attribute name="typeHL7v2" type="xs:string" use="optional" /&gt;       &lt;/xs:extension&gt;     &lt;/xs:simpleContent&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt; </pre>		

## Element t5:Unit

Namespace	http://sll-mdilab.net/T5/		
Diagram			
Type	t5:Unit		
Properties	content:	complex	
Used by	Complex Type	t5:Observation	
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>
	codingSystemName	xs:string	required
		OBX-6.3	
	id		required
		OBX-6.1	
Source	<pre> &lt;xs:element name="Unit" type="t5:Unit"&gt; </pre>		

```
</xs:element>
```

## Element t5:Timestamp

Namespace	http://sll-mdilab.net/T5/
Annotations	is extracted from parent OBR-7 or from narrowest OBX-14 (in terms of containment tree)
Diagram	<p>The diagram illustrates the relationship between the <b>t5:Timestamp</b> complex type and the <b>xs:dateTime</b> primitive type. <b>t5:Timestamp</b> is shown as a complex type that inherits from <b>xs:dateTime</b>. A note specifies that <b>t5:Timestamp</b> is extracted from parent OBR-7 or from narrowest OBX-14 (in terms of containment tree). The <b>xs:dateTime</b> type is described as a built-in primitive type representing a specific instant of time.</p>
Type	t5:Timestamp
Properties	content: complex
Used by	Complex Type t5:Observation
Source	<pre>&lt;xs:element name="Timestamp" type="t5:Timestamp"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;is extracted from parent OBR-7 or from narrowest OBX-14 (in terms of     containment tree)&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:element&gt;</pre>

## Element t5:Method

Namespace	http://sll-mdilab.net/T5/
Diagram	<p>The diagram shows the <b>t5:Method</b> element, which is a complex type.</p>
Properties	content: complex
Used by	Complex Type t5:Observation
Source	<pre>&lt;xs:element name="Method"&gt;   &lt;xs:complexType/&gt; &lt;/xs:element&gt;</pre>

## Element t5:Body\_Site

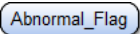
Namespace	http://sll-mdilab.net/T5/
Diagram	<p>The diagram shows the <b>t5:Body_Site</b> element, which is a complex type.</p>
Properties	content: complex
Used by	Complex Type t5:Observation
Source	<pre>&lt;xs:element name="Body_Site"&gt;   &lt;xs:complexType/&gt; &lt;/xs:element&gt;</pre>

## Element t5:Status

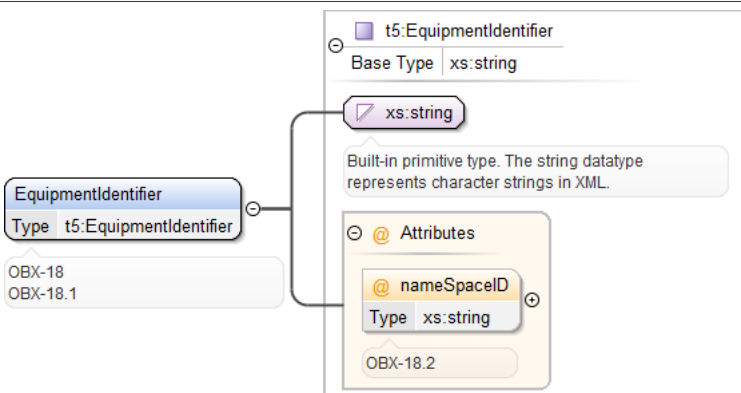
Namespace	http://sll-mdilab.net/T5/
Diagram	<p>The diagram illustrates the relationship between the <b>t5:Status</b> complex type and the <b>t5:Abnormal_Flag</b> complex type. <b>t5:Status</b> is shown as a complex type that inherits from <b>t5:Abnormal_Flag</b>.</p>
Properties	content: complex
Used by	Complex Type t5:Observation
Model	t5:Abnormal_Flag
Children	t5:Abnormal_Flag
Instance	<pre>&lt;t5:Status xmlns:t5="http://sll-mdilab.net/T5/"&gt;</pre>

	<pre>&lt;t5:Abnormal_Flag&gt;{1,1}&lt;/t5:Abnormal_Flag&gt; &lt;/t5&gt;Status&gt;</pre>
Source	<pre>&lt;xs:element name="Status"&gt;   &lt;xs:complexType&gt;     &lt;xs:sequence&gt;       &lt;xs:element ref="t5:Abnormal_Flag"/&gt;     &lt;/xs:sequence&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt;</pre>

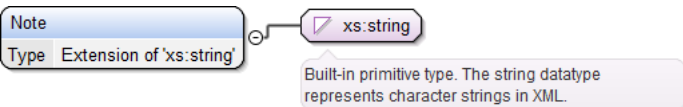
## Element t5:Abnormal\_Flag

Namespace	http://sll-mdilab.net/T5/
Diagram	
Properties	content: complex
Used by	Element t5>Status
Source	<pre>&lt;xs:element name="Abnormal_Flag"&gt;   &lt;xs:complexType/&gt; &lt;/xs:element&gt;</pre>

## Element t5:EquipmentIdentifier

Namespace	http://sll-mdilab.net/T5/														
Annotations	OBX-18 OBX-18.1														
Diagram															
Type	t5:EquipmentIdentifier														
Properties	content:	complex													
Used by	Complex Type	t5:Observation													
Attributes	<table><thead><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr></thead><tbody><tr><td>namespaceID</td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>OBX-18.2</td><td></td><td></td></tr></tbody></table>	QName	Type	Use		namespaceID	xs:string	optional			OBX-18.2				
QName	Type	Use													
namespaceID	xs:string	optional													
	OBX-18.2														
Source	<pre>&lt;xs:element name="EquipmentIdentifier" type="t5:EquipmentIdentifier"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-18&lt;/xs:documentation&gt;     &lt;xs:documentation&gt;OBX-18.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:element&gt;</pre>														

## Element t5:Note

Namespace	http://sll-mdilab.net/T5/
Diagram	

Type	extension of xs:string
Properties	content: complex
Used by	Complex Types t5:Observation, t5:Order_Observations
Source	<pre> &lt;xs:element name="Note"&gt;   &lt;xs:complexType&gt;     &lt;xs:simpleContent&gt;       &lt;xs:extension base="xs:string" /&gt;     &lt;/xs:simpleContent&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt; </pre>

## Element t5:Facet

Namespace	http://sll-mdilab.net/T5/			
Diagram				
Type	t5:Facet			
Properties	content:	complex		
Model	t5:Observation			
Children	t5:Observation			
Instance	<pre>&lt;t5:Facet id="" index="" xmlns:t5="http://sll-mdilab.net/T5/"&gt;   &lt;t5:Observation hierarchy="" index="" setid="" uid="" violatesIHE=""&gt;{1,1}&lt;/t5:Observation&gt; &lt;/t5:Facet&gt;</pre>			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:element name="Facet" type="t5:Facet"/&gt;</pre>			

## Element t5:Result

Namespace	http://sll-mdilab.net/T5/			
Diagram				
Type	extension of xs:string			
Properties	content:	complex		

Attributes	QName	Type	Use	
	<b>unit</b>		optional	
Source	<pre> &lt;xs:element name="Result"&gt;   &lt;xs:complexType&gt;     &lt;xs:simpleContent&gt;       &lt;xs:extension base="xs:string"&gt;         &lt;xs:attribute name="unit"/&gt;       &lt;/xs:extension&gt;     &lt;/xs:simpleContent&gt;   &lt;/xs:complexType&gt; &lt;/xs:element&gt; </pre>			

## Complex Type(s)

### Complex Type t5:PCD\_01\_Message

Namespace	http://sll-mdilab.net/T5/			
Annotations	HL7 v2 message as defined in IHE PCD profile <a href="http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol1.pdf">http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol1.pdf</a> <a href="http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol2.pdf">http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol2.pdf</a> <a href="http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol3.pdf">http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol3.pdf</a>			
Diagram				
Used by	Element	t5:PCD_01_Message		
Model	t5:Sending_Application , t5:Sending_Facility , t5:Patient_Result+			
Children	t5:Patient_Result, t5:Sending_Application, t5:Sending_Facility			
Attributes	QName	Type	Use	
	id	xs:string	optional	
		MSH-10		
	timeStam	xs:dateTime	optional	
		MSH-7		
Source	<pre>&lt;xs:complexType name="PCD_01_Message"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;HL7 v2 message as defined in IHE PCD profile <a href="http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol1.pdf">http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol1.pdf</a> <a href="http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol2.pdf">http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol2.pdf</a> <a href="http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol3.pdf">http://ihe.net/uploadedFiles/Documents/PCD/IHE_PCD_TF_Vol3.pdf</a>&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:sequence&gt;     &lt;xs:element name="Sending_Application"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;MSH-3 MSH-3.1 is used as Application ID, which is an identifier of the System&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;       &lt;xs:complexType&gt;         &lt;xs:attribute name="idLocal" use="optional"&gt;</pre>			

```

        <xs:annotation>
          <xs:documentation>MSH-3.1</xs:documentation>
        </xs:annotation>
      </xs:attribute>
      <xs:attribute name="idUniversal">
        <xs:annotation>
          <xs:documentation>MSH-3.2</xs:documentation>
        </xs:annotation>
      </xs:attribute>
      <xs:attribute name="idUniversalType">
        <xs:annotation>
          <xs:documentation>MSH-3.3</xs:documentation>
        </xs:annotation>
      </xs:attribute>
    </xs:complexType>
  </xs:element>
  <xs:element name="Sending_Facility">
    <xs:annotation>
      <xs:documentation>MSH-4</xs:documentation>
      <xs:documentation>to distinguish Hospital/department</xs:documentation>
    </xs:annotation>
    <xs:complexType>
      <xs:attribute name="idLocal" use="optional">
        <xs:annotation>
          <xs:documentation>MSH-4.1</xs:documentation>
        </xs:annotation>
      </xs:attribute>
      <xs:attribute name="idUniversal">
        <xs:annotation>
          <xs:documentation>MSH-4.2</xs:documentation>
        </xs:annotation>
      </xs:attribute>
      <xs:attribute name="idUniversalType">
        <xs:annotation>
          <xs:documentation>MSH-4.3</xs:documentation>
        </xs:annotation>
      </xs:attribute>
    </xs:complexType>
  </xs:element>
  <xs:element ref="t5:Patient_Result" maxOccurs="unbounded" />
</xs:sequence>
<xs:attribute name="timeStamp" type="xs:dateTime">
  <xs:annotation>
    <xs:documentation>MSH-7</xs:documentation>
  </xs:annotation>
</xs:attribute>
<xs:attribute name="id" type="xs:string">
  <xs:annotation>
    <xs:documentation>MSH-10</xs:documentation>
  </xs:annotation>
</xs:attribute>
</xs:complexType>

```

## Complex Type t5:Patient\_Result

Namespace	http://sll-mdilab.net/T5/
Diagram	
Used by	Element t5:Patient_Result
Model	t5:Patient , t5:Order_Observations+
Children	t5:Order_Observations, t5:Patient
Source	<pre> &lt;xs:complexType name="Patient_Result"&gt;   &lt;xs:sequence&gt;     &lt;xs:element ref="t5:Patient"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;PID&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Order_Observations" maxOccurs="unbounded" /&gt;   &lt;/xs:sequence&gt; &lt;/xs:complexType&gt; </pre>

**Complex Type t5:Patient**

Namespace	http://sll-mdilab.net/T5/		
Diagram	<p>The diagram illustrates the structure of the <code>t5:Patient</code> complex type. It is a sequence of three elements: <code>t5:Identifier</code> (1..6), <code>Name</code>, and <code>t5:Location</code> (0..1). <code>t5:Identifier</code> is an extension of <code>xs:string</code> and has a documentation of <code>PID-3</code>. <code>Name</code> has a documentation of <code>PID-5</code>. <code>t5:Location</code> is of type <code>t5:Location</code> and has a documentation of <code>PV1-3</code>.</p>		
Used by	Element	<code>t5:Patient</code>	
Model	<code>t5:Identifier{1,6} , t5:Name , t5:Location{0,1}</code>		
Children	<code>t5:Identifier, t5:Location, t5:Name</code>		
Source	<pre>&lt;xs:complexType name="Patient"&gt;   &lt;xs:sequence minOccurs="0"&gt;     &lt;xs:element ref="t5:Identifier" maxOccurs="6"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;PID-3&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element name="Name"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;PID-5&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element minOccurs="0" ref="t5:Location"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;PV1-3&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;   &lt;/xs:sequence&gt; &lt;/xs:complexType&gt;</pre>		

**Complex Type t5:Location**

Namespace	http://sll-mdilab.net/T5/		
Diagram			
Used by	Element	t5:Location	
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>
	<b>bed</b>	xs:string	optional
		PV1-3.3	
	<b>pointOfCare</b>	xs:string	optional
		PV1-3.1	
	<b>room</b>	xs:string	optional
		PV1-3.2	



Source	<pre> &lt;xs:complexType name="Location"&gt;   &lt;xs:attribute name="room" type="xs:string"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;PV1-3.2&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt;   &lt;xs:attribute name="bed" use="optional" type="xs:string"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;PV1-3.3&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt;   &lt;xs:attribute name="pointOfCare" type="xs:string" use="optional"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;PV1-3.1&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt; </pre>
--------	---

## Complex Type t5:Order\_Observations

Namespace	http://sll-mdilab.net/T5/
Diagram	
Used by	Element t5:Order_Observations
Model	t5:Order{0,1} , t5:MDS* , t5:Note*
Children	t5:MDS, t5:Note, t5:Order
Source	<pre> &lt;xs:complexType name="Order_Observations"&gt;   &lt;xs:sequence minOccurs="0"&gt;     &lt;xs:element ref="t5:Order" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:MDS" maxOccurs="unbounded" minOccurs="0"/&gt;     &lt;xs:element ref="t5:Note" maxOccurs="unbounded" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;NTE&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;   &lt;/xs:sequence&gt; &lt;/xs:complexType&gt; </pre>

## Complex Type t5:Order

Namespace	http://sll-mdilab.net/T5/
Diagram	

Used by	Element t5:Order			
Model	t5:FillerOrderNumber , t5:UniversalServiceID			
Children	t5:FillerOrderNumber, t5:UniversalServiceID			
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>	
	<b>timeStamp</b>	xs:dateTime	optional	
		OBR-7		
Source	<pre> &lt;xs:complexType name="Order"&gt;   &lt;xs:sequence&gt;     &lt;xs:element name="FillerOrderNumber" maxOccurs="1"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-3.1&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;       &lt;xs:complexType&gt;         &lt;xs:attribute name="idLocal" use="optional"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;OBR-3.1&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;         &lt;xs:attribute name="nameSpaceID" type="xs:string"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;OBR-3.2&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;         &lt;xs:attribute name="idUniversal"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;OBR-3.3&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;         &lt;xs:attribute name="idUniversalType"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;MSH-3.4&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;       &lt;/xs:complexType&gt;     &lt;/xs:element&gt;     &lt;xs:element name="UniversalServiceID" maxOccurs="1" minOccurs="1"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBR-4.2 text&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;       &lt;xs:complexType&gt;         &lt;xs:attribute name="id" use="optional"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;OBR-4.1 id&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;         &lt;xs:attribute name="codingSystemName" type="xs:string" use="optional"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;OBR-4.3&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;       &lt;/xs:complexType&gt;     &lt;/xs:element&gt;   &lt;/xs:sequence&gt;   &lt;xs:attribute name="timeStamp" type="xs:dateTime"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;OBR-7&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt; </pre>			

## Complex Type t5:MDS

Namespace	http://sll-mdilab.net/T5/
-----------	---------------------------

Diagram				
Used by	Element	t5:MDS		
Model	t5:VMD{0,1} , t5:Observation*			
Children	t5:Observation, t5:VMD			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.0.0.0		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:complexType name="MDS"&gt;   &lt;xs:sequence maxOccurs="unbounded"&gt;     &lt;xs:element ref="t5:VMD" maxOccurs="1" minOccurs="0"/&gt;     &lt;xs:element ref="t5:Observation" maxOccurs="unbounded" minOccurs="0"/&gt;   &lt;/xs:sequence&gt;   &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt;   &lt;xs:attribute name="id"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.0.0.0&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt;</pre>			

## Complex Type t5:VMD

Namespace	http://sll-mdilab.net/T5/		
Diagram			
Used by	Element	t5:VMD	
Model	t5:CHAN{0,1} , t5:Observation*		

Children	t5:CHAN, t5:Observation			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.y.0.0		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:complexType name="VMD"&gt;   &lt;xs:sequence maxOccurs="unbounded"&gt;     &lt;xs:element ref="t5:CHAN" minOccurs="0"/&gt;     &lt;xs:element ref="t5:Observation" maxOccurs="unbounded" minOccurs="0"/&gt;   &lt;/xs:sequence&gt;   &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt;   &lt;xs:attribute name="id"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.0.0&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt;</pre>			

## Complex Type t5:CHAN

Namespace	http://sll-mdilab.net/T5/																				
Diagram																					
Used by	Element	t5:CHAN																			
Model	t5:Metric{0,1} , t5:Observation*																				
Children	t5:Metric, t5:Observation																				
Attributes	<table><thead><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr></thead><tbody><tr><td>id</td><td></td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">OBX-3.2 extracted from segment where OBX-4 is x.y.z.0</td></tr><tr><td>index</td><td>xs:int</td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">index in ContainmentTree, extracted from OBX-4</td></tr></tbody></table>	QName	Type	Use		id		optional			OBX-3.2 extracted from segment where OBX-4 is x.y.z.0			index	xs:int	optional			index in ContainmentTree, extracted from OBX-4		
QName	Type	Use																			
id		optional																			
	OBX-3.2 extracted from segment where OBX-4 is x.y.z.0																				
index	xs:int	optional																			
	index in ContainmentTree, extracted from OBX-4																				
Source	<pre>&lt;xs:complexType name="CHAN"&gt;   &lt;xs:sequence maxOccurs="unbounded"&gt;     &lt;xs:element ref="t5:Metric" minOccurs="0"/&gt;     &lt;xs:element ref="t5:Observation" maxOccurs="unbounded" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;   &lt;/xs:sequence&gt;   &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt;</pre>																				

```

</xs:attribute>
<xs:attribute name="id">
  <xs:annotation>
    <xs:documentation>OBX-3.2 extracted from segment where OBX-4 is x.y.z.0</xs:documentation>
  </xs:annotation>
</xs:attribute>
</xs:complexType>

```

## Complex Type t5:Metric

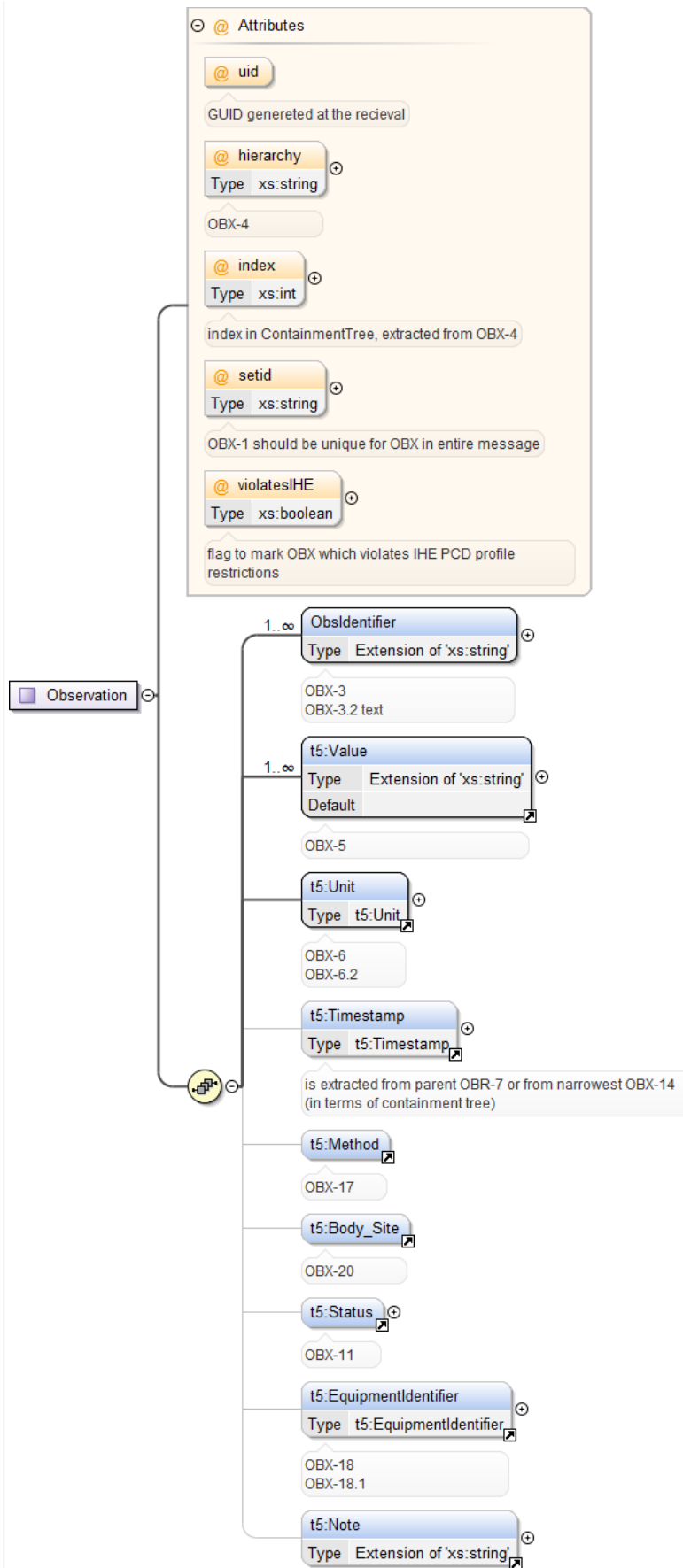
Namespace	http://sll-mdilab.net/T5/			
Diagram				
Used by	Element	t5:Metric		
Model	t5:Facet{0,1} , t5:Observation*			
Children	t5:Facet, t5:Observation			
Attributes	QName	Type	Use	
	id		optional	
		OBX-3.2 extracted from segment where OBX-4 is x.y.z.w		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
Source	<pre>&lt;xs:complexType name="Metric"&gt;   &lt;xs:sequence maxOccurs="unbounded"&gt;     &lt;xs:element name="Facet" minOccurs="0"&gt;       &lt;xs:complexType&gt;         &lt;xs:sequence&gt;           &lt;xs:element ref="t5:Observation"/&gt;         &lt;/xs:sequence&gt;         &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;         &lt;xs:attribute name="id"&gt;           &lt;xs:annotation&gt;             &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t&lt;/ xs:documentation&gt;           &lt;/xs:annotation&gt;         &lt;/xs:attribute&gt;       &lt;/xs:complexType&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Observation" maxOccurs="unbounded" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;   &lt;/xs:sequence&gt;   &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt;   &lt;xs:attribute name="id"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt;</pre>			

	<pre>&lt;/xs:annotation&gt; &lt;/xs:attribute&gt; &lt;/xs:complexType&gt;</pre>
--	---

### Complex Type t5:Observation

Namespace	<a href="http://sll-mdilab.net/T5/">http://sll-mdilab.net/T5/</a>
-----------	---

Diagram



Used by

Element

t5:Observation

Model	t5:ObsIdentifier+, t5:Value+, t5:Unit, t5:Timestamp{0,1}, t5:Method{0,1}, t5:Body_Site{0,1}, t5:Status{0,1}, t5:EquipmentIdentifier{0,1}, t5:Note{0,1}			
Children	t5:Body_Site, t5:EquipmentIdentifier, t5:Method, t5:Note, t5:ObsIdentifier, t5:Status, t5:Timestamp, t5:Unit, t5:Value			
Attributes	QName	Type	Use	
	hierarchy	xs:string	optional	
		OBX-4		
	index	xs:int	optional	
		index in ContainmentTree, extracted from OBX-4		
	setid	xs:string	optional	
		OBX-1 should be unique for OBX in entire message		
	uid		optional	
		GUID genereted at the recieval		
	violatesIHE	xs:boolean	optional	
	flag to mark OBX which violates IHE PCD profile restrictions			
Source	<pre>&lt;xs:complexType name="Observation"&gt;   &lt;xs:sequence&gt;     &lt;xs:element name="ObsIdentifier" maxOccurs="unbounded"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX-3&lt;/xs:documentation&gt;         &lt;xs:documentation&gt;OBX-3.2 text&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;       &lt;xs:complexType&gt;         &lt;xs:simpleContent&gt;           &lt;xs:extension base="xs:string"&gt;             &lt;xs:attribute name="id" use="required"&gt;               &lt;xs:annotation&gt;                 &lt;xs:documentation&gt;OBX-3.1&lt;/xs:documentation&gt;               &lt;/xs:annotation&gt;             &lt;/xs:attribute&gt;             &lt;xs:attribute name="codingSystemName" type="xs:string" use="required"&gt;               &lt;xs:annotation&gt;                 &lt;xs:documentation&gt;OBX-3.3&lt;/xs:documentation&gt;               &lt;/xs:annotation&gt;             &lt;/xs:attribute&gt;             &lt;xs:attribute name="isAlternate" type="xs:boolean"/&gt;             &lt;xs:attribute name="notFound" type="xs:boolean"&gt;               &lt;xs:annotation&gt;                 &lt;xs:documentation&gt;flag to be set if the identifier term could not be found in the nomenclature (i.e. can not be found in IEEE/IHE Rosetta table)&lt;/xs:documentation&gt;               &lt;/xs:annotation&gt;             &lt;/xs:attribute&gt;           &lt;/xs:extension&gt;         &lt;/xs:simpleContent&gt;       &lt;/xs:complexType&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Value" maxOccurs="unbounded"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX-5&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Unit"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX-6&lt;/xs:documentation&gt;         &lt;xs:documentation&gt;OBX-6.2&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Timestamp" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;is extracted from parent OBR-7 or from narrowest OBX-14 (in terms of containment tree)&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Method" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX-17&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Body_Site" minOccurs="0"&gt;       &lt;xs:annotation&gt;         &lt;xs:documentation&gt;OBX-20&lt;/xs:documentation&gt;       &lt;/xs:annotation&gt;     &lt;/xs:element&gt;     &lt;xs:element ref="t5:Status" minOccurs="0"&gt;       &lt;xs:annotation&gt;</pre>			



```

        <xs:documentation>OBX-11</xs:documentation>
      </xs:annotation>
    </xs:element>
    <xs:element minOccurs="0" ref="t5:EquipmentIdentifier"/>
    <xs:element ref="t5:Note" minOccurs="0"/>
  </xs:sequence>
  <xs:attribute name="uid">
    <xs:annotation>
      <xs:documentation>GUID generated at the recieval</xs:documentation>
    </xs:annotation>
  </xs:attribute>
  <xs:attribute name="hierarchy" type="xs:string">
    <xs:annotation>
      <xs:documentation>OBX-4</xs:documentation>
    </xs:annotation>
  </xs:attribute>
  <xs:attribute name="index" form="unqualified" type="xs:int">
    <xs:annotation>
      <xs:documentation>index in ContainmentTree, extracted from OBX-4</xs:documentation>
    </xs:annotation>
  </xs:attribute>
  <xs:attribute name="setid" type="xs:string">
    <xs:annotation>
      <xs:documentation>OBX-1 should be unique for OBX in entire message</xs:documentation>
    </xs:annotation>
  </xs:attribute>
  <xs:attribute name="violatesIHE" type="xs:boolean" use="optional">
    <xs:annotation>
      <xs:documentation>flag to mark OBX which violates IHE PCD profile restrictions</
xs:documentation>
    </xs:annotation>
  </xs:attribute>
</xs:complexType>

```

## Complex Type t5:Unit

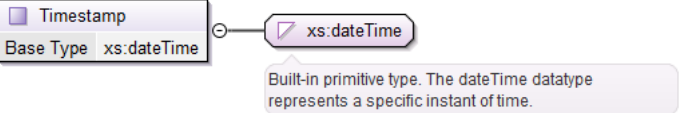
Namespace	http://sll-mdilab.net/T5/		
Diagram	<p>The diagram illustrates the complex type <code>t5:Unit</code> as an extension of the <code>xs:string</code> primitive type. A box labeled <code>Unit</code> with <code>Base Type xs:string</code> is connected to the <code>xs:string</code> type. The <code>Unit</code> type has two attributes: <code>@ id</code> (required, OBX-6.1) and <code>@ codingSystemName</code> (required, OBX-6.3, of type <code>xs:string</code>). A callout for <code>xs:string</code> states: "Built-in primitive type. The string datatype represents character strings in XML."</p>		
Type	extension of xs:string		
Used by	Element	t5:Unit	
Attributes	<b>QName</b>	<b>Type</b>	<b>Use</b>
	<b>codingSystemName</b>	xs:string	required
		OBX-6.3	
	<b>id</b>		required
		OBX-6.1	
Source	<pre> &lt;xs:complexType name="Unit"&gt;   &lt;xs:simpleContent&gt;     &lt;xs:extension base="xs:string"&gt;       &lt;xs:attribute name="id" use="required"&gt;         &lt;xs:annotation&gt;           &lt;xs:documentation&gt;OBX-6.1&lt;/xs:documentation&gt;         &lt;/xs:annotation&gt;       &lt;/xs:attribute&gt;       &lt;xs:attribute name="codingSystemName" type="xs:string" use="required"&gt;         &lt;xs:annotation&gt;           &lt;xs:documentation&gt;OBX-6.3&lt;/xs:documentation&gt;         &lt;/xs:annotation&gt;       &lt;/xs:attribute&gt;     &lt;/xs:extension&gt;   &lt;/xs:simpleContent&gt; &lt;/xs:complexType&gt; </pre>		

```

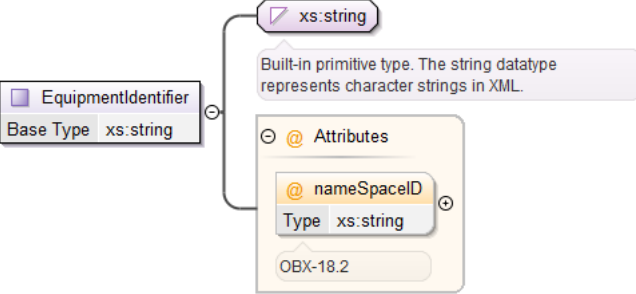
</xs:extension>
</xs:simpleContent>
</xs:complexType>

```

### Complex Type t5:Timestamp

Namespace	http://sll-mdilab.net/T5/
Diagram	
Type	extension of xs:dateTime
Used by	Element t5:Timestamp
Source	<pre> &lt;xs:complexType name="Timestamp"&gt;   &lt;xs:simpleContent&gt;     &lt;xs:extension base="xs:dateTime" /&gt;   &lt;/xs:simpleContent&gt; &lt;/xs:complexType&gt; </pre>

### Complex Type t5:EquipmentIdentifier

Namespace	http://sll-mdilab.net/T5/														
Diagram															
Type	extension of xs:string														
Used by	Element t5:EquipmentIdentifier														
Attributes	<table><thead><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr></thead><tbody><tr><td>namespaceID</td><td>xs:string</td><td>optional</td><td></td></tr><tr><td></td><td>OBX-18.2</td><td></td><td></td></tr></tbody></table>	QName	Type	Use		namespaceID	xs:string	optional			OBX-18.2				
QName	Type	Use													
namespaceID	xs:string	optional													
	OBX-18.2														
Source	<pre>&lt;xs:complexType name="EquipmentIdentifier"&gt;   &lt;xs:simpleContent&gt;     &lt;xs:extension base="xs:string"&gt;       &lt;xs:attribute name="namespaceID" type="xs:string"&gt;         &lt;xs:annotation&gt;           &lt;xs:documentation&gt;OBX-18.2&lt;/xs:documentation&gt;         &lt;/xs:annotation&gt;       &lt;/xs:attribute&gt;     &lt;/xs:extension&gt;   &lt;/xs:simpleContent&gt; &lt;/xs:complexType&gt;</pre>														

### Complex Type t5:Facet

Namespace	http://sll-mdilab.net/T5/
Annotations	This Device Containment Tree component/level is defined in IHE PCD WCM supplement <a href="http://ihe.net/Technical_Framework/upload/IHE_PCD_Suppl_WCM.pdf">http://ihe.net/Technical_Framework/upload/IHE_PCD_Suppl_WCM.pdf</a>

Diagram	<p>The diagram illustrates the structure of the <code>t5:Facet</code> complex type. It features a main container box labeled <code>@ Attributes</code> which holds two attributes: <code>@ index</code> (of type <code>xs:int</code>) and <code>@ id</code>. The <code>@ index</code> attribute is annotated with the text "index in ContainmentTree, extracted from OBX-4". The <code>@ id</code> attribute is annotated with "OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t". A separate box labeled <code>Facet</code> is connected to the main container and contains the text "This Device Containment Tree component/level is defined in IHE PCD WCM supplement...". The <code>Facet</code> box is also connected to a box representing the <code>t5:Observation</code> element, which has a type of <code>t5:Observation</code>.</p>																				
Used by	Element	t5:Facet																			
Model	t5:Observation																				
Children	t5:Observation																				
Attributes	<table><thead><tr><th>QName</th><th>Type</th><th>Use</th><th></th></tr></thead><tbody><tr><td>id</td><td></td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t</td></tr><tr><td>index</td><td>xs:int</td><td>optional</td><td></td></tr><tr><td></td><td colspan="3">index in ContainmentTree, extracted from OBX-4</td></tr></tbody></table>	QName	Type	Use		id		optional			OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t			index	xs:int	optional			index in ContainmentTree, extracted from OBX-4		
QName	Type	Use																			
id		optional																			
	OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t																				
index	xs:int	optional																			
	index in ContainmentTree, extracted from OBX-4																				
Source	<pre>&lt;xs:complexType name="Facet"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;This Device Containment Tree component/level is defined in IHE PCD WCM     supplement http://ihe.net/Technical_Framework/upload/IHE_PCD_Suppl_WCM.pdf&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;   &lt;xs:sequence&gt;     &lt;xs:element ref="t5:Observation"/&gt;   &lt;/xs:sequence&gt;   &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt;   &lt;xs:attribute name="id"&gt;     &lt;xs:annotation&gt;       &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t&lt;/xs:documentation&gt;     &lt;/xs:annotation&gt;   &lt;/xs:attribute&gt; &lt;/xs:complexType&gt;</pre>																				

**Namespace:** ""

### Attribute(s)

**Attribute t5:PCD\_01\_Message / t5:Sending\_Application / @idLocal**

Namespace	No namespace
Annotations	MSH-3.1
Properties	use: optional
Used by	Element t5:PCD_01_Message/t5:Sending_Application
Source	<pre> &lt;xs:attribute name="idLocal" use="optional"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-3.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:PCD\_01\_Message / t5:Sending\_Application / @idUniversal**

Namespace	No namespace
Annotations	MSH-3.2

Used by	Element t5:PCD_01_Message/t5:Sending_Application
Source	<pre>&lt;xs:attribute name="idUniversal"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-3.2&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:PCD\_01\_Message / t5:Sending\_Application / @idUniversalType**

Namespace	No namespace
Annotations	MSH-3.3
Used by	Element t5:PCD_01_Message/t5:Sending_Application
Source	<pre>&lt;xs:attribute name="idUniversalType"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-3.3&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:PCD\_01\_Message / t5:Sending\_Facility / @idLocal**

Namespace	No namespace
Annotations	MSH-4.1
Properties	use: optional
Used by	Element t5:PCD_01_Message/t5:Sending_Facility
Source	<pre>&lt;xs:attribute name="idLocal" use="optional"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-4.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:PCD\_01\_Message / t5:Sending\_Facility / @idUniversal**

Namespace	No namespace
Annotations	MSH-4.2
Used by	Element t5:PCD_01_Message/t5:Sending_Facility
Source	<pre>&lt;xs:attribute name="idUniversal"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-4.2&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:PCD\_01\_Message / t5:Sending\_Facility / @idUniversalType**

Namespace	No namespace
Annotations	MSH-4.3
Used by	Element t5:PCD_01_Message/t5:Sending_Facility
Source	<pre>&lt;xs:attribute name="idUniversalType"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-4.3&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Identifier / @authority**

Namespace	No namespace
Annotations	PID-3.4
Properties	use: required
Used by	Element t5:Identifier
Source	<pre>&lt;xs:attribute name="authority" use="required"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PID-3.4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt;</pre>

	</xs:attribute>
--	-----------------

**Attribute t5:Identifier / @typeCode**

Namespace	No namespace
Annotations	PID-3.5
Used by	Element t5:Identifier
Source	<pre>&lt;xs:attribute name="typeCode"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PID-3.5&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Location / @room**

Namespace	No namespace
Annotations	PV1-3.2
Type	xs:string
Properties	content: simple
Used by	Complex Type t5:Location
Source	<pre>&lt;xs:attribute name="room" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PV1-3.2&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Location / @bed**

Namespace	No namespace
Annotations	PV1-3.3
Type	xs:string
Properties	use: optional
Used by	Complex Type t5:Location
Source	<pre>&lt;xs:attribute name="bed" use="optional" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PV1-3.3&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Location / @pointOfCare**

Namespace	No namespace
Annotations	PV1-3.1
Type	xs:string
Properties	use: optional
Used by	Complex Type t5:Location
Source	<pre>&lt;xs:attribute name="pointOfCare" type="xs:string" use="optional"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;PV1-3.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Order / t5:FillerOrderNumber / @idLocal**

Namespace	No namespace
Annotations	OBR-3.1
Properties	use: optional
Used by	Element t5:Order/t5:FillerOrderNumber
Source	<pre>&lt;xs:attribute name="idLocal" use="optional"&gt;</pre>

```

<xs:annotation>
  <xs:documentation>OBR-3.1</xs:documentation>
</xs:annotation>
</xs:attribute>

```

**Attribute t5:Order / t5:FillerOrderNumber / @nameSpaceID**

Namespace	No namespace
Annotations	OBR-3.2
Type	xs:string
Properties	content: simple
Used by	Element t5:Order/t5:FillerOrderNumber
Source	<pre> &lt;xs:attribute name="nameSpaceID" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBR-3.2&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Order / t5:FillerOrderNumber / @idUniversal**

Namespace	No namespace
Annotations	OBR-3.3
Used by	Element t5:Order/t5:FillerOrderNumber
Source	<pre> &lt;xs:attribute name="idUniversal"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBR-3.3&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Order / t5:FillerOrderNumber / @idUniversalType**

Namespace	No namespace
Annotations	MSH-3.4
Used by	Element t5:Order/t5:FillerOrderNumber
Source	<pre> &lt;xs:attribute name="idUniversalType"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-3.4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Order / t5:UniversalServiceID / @id**

Namespace	No namespace
Annotations	OBR-4.1 id
Properties	use: optional
Used by	Element t5:Order/t5:UniversalServiceID
Source	<pre> &lt;xs:attribute name="id" use="optional"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBR-4.1 id&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Order / t5:UniversalServiceID / @codingSystemName**

Namespace	No namespace
Annotations	OBR-4.3
Type	xs:string
Properties	use: optional
Used by	Element t5:Order/t5:UniversalServiceID
Source	<pre> &lt;xs:attribute name="codingSystemName" type="xs:string" use="optional"&gt;   &lt;xs:annotation&gt; </pre>

```

<xs:documentation>OBR-4.3</xs:documentation>
</xs:annotation>
</xs:attribute>

```

**Attribute t5:Order / @timeStamp**

Namespace	No namespace
Annotations	OBR-7
Type	xs:dateTime
Properties	content: simple
Used by	Complex Type t5:Order
Source	<pre> &lt;xs:attribute name="timeStamp" type="xs:dateTime"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBR-7&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Observation / t5:ObsIdentifier / @id**

Namespace	No namespace
Annotations	OBX-3.1
Properties	use: required
Used by	Element t5:Observation/t5:ObsIdentifier
Source	<pre> &lt;xs:attribute name="id" use="required"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Observation / t5:ObsIdentifier / @codingSystemName**

Namespace	No namespace
Annotations	OBX-3.3
Type	xs:string
Properties	use: required
Used by	Element t5:Observation/t5:ObsIdentifier
Source	<pre> &lt;xs:attribute name="codingSystemName" type="xs:string" use="required"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.3&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Observation / t5:ObsIdentifier / @isAlternate**

Namespace	No namespace
Type	xs:boolean
Properties	content: simple
Used by	Element t5:Observation/t5:ObsIdentifier
Source	<pre> &lt;xs:attribute name="isAlternate" type="xs:boolean"/&gt; </pre>

**Attribute t5:Observation / t5:ObsIdentifier / @notFound**

Namespace	No namespace
Annotations	flag to be set if the identifier term could not be found in the nomenclature (i.e. can not be found in IEEE/IHE Rosetta table)
Type	xs:boolean
Properties	content: simple
Used by	Element t5:Observation/t5:ObsIdentifier

Source	<pre>&lt;xs:attribute name="notFound" type="xs:boolean"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;flag to be set if the identifier term could not be found in the nomenclature     (i.e. can not be found in IEEE/IHE Rosetta table)&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>
--------	--

**Attribute t5:Value / @typeHL7v2**

Namespace	No namespace	
Type	xs:string	
Properties	use:	optional
Used by	Element	t5:Value
Source	<pre>&lt;xs:attribute name="typeHL7v2" type="xs:string" use="optional"/&gt;</pre>	

**Attribute t5:Unit / @id**

Namespace	No namespace	
Annotations	OBX-6.1	
Properties	use:	required
Used by	Complex Type	t5:Unit
Source	<pre>&lt;xs:attribute name="id" use="required"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-6.1&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>	

**Attribute t5:Unit / @codingSystemName**

Namespace	No namespace	
Annotations	OBX-6.3	
Type	xs:string	
Properties	use:	required
Used by	Complex Type	t5:Unit
Source	<pre>&lt;xs:attribute name="codingSystemName" type="xs:string" use="required"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-6.3&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>	

**Attribute t5:EquipmentIdentifier / @nameSpaceID**

Namespace	No namespace	
Annotations	OBX-18.2	
Type	xs:string	
Properties	content:	simple
Used by	Complex Type	t5:EquipmentIdentifier
Source	<pre>&lt;xs:attribute name="nameSpaceID" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-18.2&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>	

**Attribute t5:Observation / @uid**

Namespace	No namespace	
Annotations	GUID genereted at the recieval	
Used by	Complex Type	t5:Observation
Source	<pre>&lt;xs:attribute name="uid"&gt;</pre>	



```

<xs:annotation>
  <xs:documentation>GUID generated at the recieval</xs:documentation>
</xs:annotation>
</xs:attribute>

```

**Attribute t5:Observation / @hierarchy**

Namespace	No namespace
Annotations	OBX-4
Type	xs:string
Properties	content: simple
Used by	Complex Type t5:Observation
Source	<pre> &lt;xs:attribute name="hierarchy" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Observation / @index**

Namespace	No namespace
Annotations	index in ContainmentTree, extracted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Complex Type t5:Observation
Source	<pre> &lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Observation / @setid**

Namespace	No namespace
Annotations	OBX-1 should be unique for OBX in entire message
Type	xs:string
Properties	content: simple
Used by	Complex Type t5:Observation
Source	<pre> &lt;xs:attribute name="setid" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-1 should be unique for OBX in entire message&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Observation / @violatesIHE**

Namespace	No namespace
Annotations	flag to mark OBX which violates IHE PCD profile restrictions
Type	xs:boolean
Properties	use: optional
Used by	Complex Type t5:Observation
Source	<pre> &lt;xs:attribute name="violatesIHE" type="xs:boolean" use="optional"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;flag to mark OBX which violates IHE PCD profile restrictions&lt;/ xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt; </pre>

**Attribute t5:Metric / t5:Facet / @index**

Namespace	No namespace
-----------	--------------

Annotations	index in ContainmentTree, extracted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Element t5:Metric/t5:Facet
Source	<pre>&lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Metric / t5:Facet / @id**

Namespace	No namespace
Annotations	OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t
Used by	Element t5:Metric/t5:Facet
Source	<pre>&lt;xs:attribute name="id"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Metric / @index**

Namespace	No namespace
Annotations	index in ContainmentTree, extracted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Complex Type t5:Metric
Source	<pre>&lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Metric / @id**

Namespace	No namespace
Annotations	OBX-3.2 extracted from segment where OBX-4 is x.y.z.w
Used by	Complex Type t5:Metric
Source	<pre>&lt;xs:attribute name="id"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:CHAN / @index**

Namespace	No namespace
Annotations	index in ContainmentTree, extracted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Complex Type t5:CHAN
Source	<pre>&lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:CHAN / @id**

Namespace	No namespace
-----------	--------------

Annotations	OBX-3.2 extracted from segment where OBX-4 is x.y.z.0
Used by	Complex Type t5:CHAN
Source	<pre>&lt;xs:attribute name="id"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.0&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:VMD / @index**

Namespace	No namespace
Annotations	index in ContainmentTree, extarcted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Complex Type t5:VMD
Source	<pre>&lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extarcted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:VMD / @id**

Namespace	No namespace
Annotations	OBX-3.2 extracted from segment where OBX-4 is x.y.0.0
Used by	Complex Type t5:VMD
Source	<pre>&lt;xs:attribute name="id"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.0.0&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:MDS / @index**

Namespace	No namespace
Annotations	index in ContainmentTree, extarcted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Complex Type t5:MDS
Source	<pre>&lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extarcted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:MDS / @id**

Namespace	No namespace
Annotations	OBX-3.2 extracted from segment where OBX-4 is x.0.0.0
Used by	Complex Type t5:MDS
Source	<pre>&lt;xs:attribute name="id"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.0.0.0&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:PCD\_01\_Message / @timeStamp**

Namespace	No namespace
Annotations	MSH-7
Type	xs:dateTime

Properties	content: simple
Used by	Complex Type t5:PCD_01_Message
Source	<pre>&lt;xs:attribute name="timeStamp" type="xs:dateTime"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-7&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:PCD\_01\_Message / @id**

Namespace	No namespace
Annotations	MSH-10
Type	xs:string
Properties	content: simple
Used by	Complex Type t5:PCD_01_Message
Source	<pre>&lt;xs:attribute name="id" type="xs:string"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;MSH-10&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Facet / @index**

Namespace	No namespace
Annotations	index in ContainmentTree, extracted from OBX-4
Type	xs:int
Properties	content: simple
Used by	Complex Type t5:Facet
Source	<pre>&lt;xs:attribute name="index" form="unqualified" type="xs:int"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;index in ContainmentTree, extracted from OBX-4&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Facet / @id**

Namespace	No namespace
Annotations	OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t
Used by	Complex Type t5:Facet
Source	<pre>&lt;xs:attribute name="id"&gt;   &lt;xs:annotation&gt;     &lt;xs:documentation&gt;OBX-3.2 extracted from segment where OBX-4 is x.y.z.w.t&lt;/xs:documentation&gt;   &lt;/xs:annotation&gt; &lt;/xs:attribute&gt;</pre>

**Attribute t5:Result / @unit**

Namespace	No namespace
Used by	Element t5:Result
Source	<pre>&lt;xs:attribute name="unit"/&gt;</pre>